


In Vitro Diagnostic Reagent

Hemolysing Agent–MNR

Used prior to blood cell analysis to lyse red blood cells, release hemoglobin and preserve the morphology of the cells to be analyzed, thus facilitating cell differentiation and counting or quantitative determination of hemoglobin

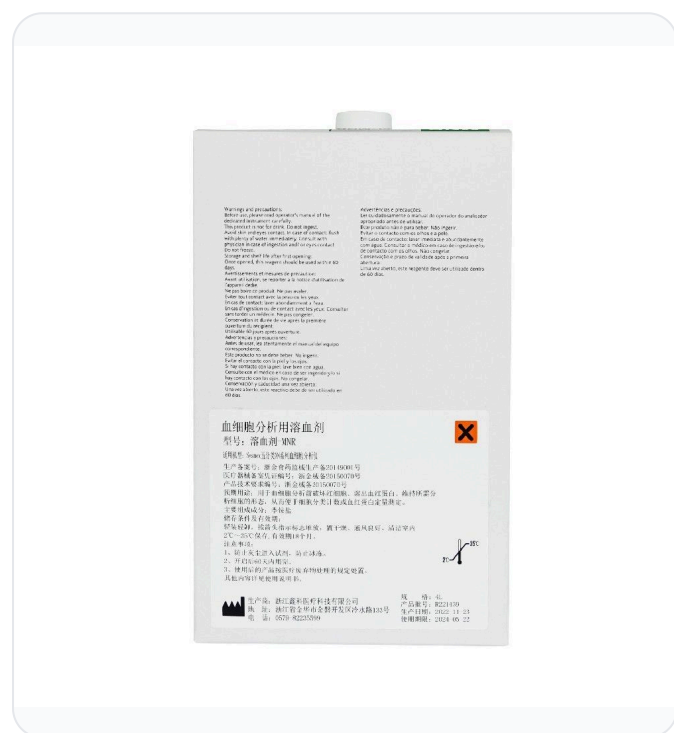

Product Photo

Product Model

Hemolysing Agent–MNR

Packaging Specifications

4 L, 5 L, 20 L

Product Code

p1095

Intended Use

Used prior to blood cell analysis to lyse red blood cells, release hemoglobin and preserve the morphology of the cells to be analyzed, thus facilitating cell differentiation and counting or quantitative determination of hemoglobin

Principle

It lyses red blood cells while keeping the morphology of white blood cells unchanged, enabling the analyzer to count and classify white blood cells by the laser light scatter method

Composition

Quaternary ammonium salt...

