


## In Vitro Diagnostic Reagent

## Hemolysing Agent-86D

Used prior to blood cell analysis to lyse red blood cells, release hemoglobin and preserve the morphology of the cells to be analyzed; when required, the white blood cells are stained before cell differentiation and counting, or cell differentiation and counting and/or quantitative determination of hemoglobin are carried out directly

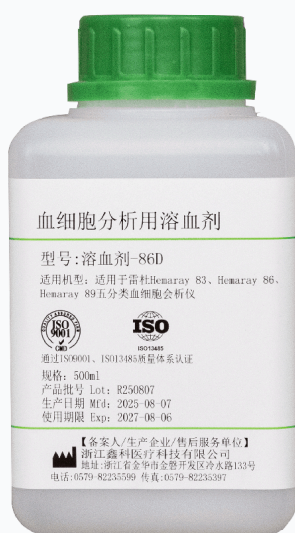

Product Photo

|                          |                         |
|--------------------------|-------------------------|
| Product Model            | Hemolysing Agent-86D    |
| Packaging Specifications | 500 mL, 5 L, 10 L, 20 L |
| Product Code             | p1140                   |

## Intended Use

Used prior to blood cell analysis to lyse red blood cells, release hemoglobin and preserve the morphology of the cells to be analyzed; when required, the white blood cells are stained before cell differentiation and counting, or cell differentiation and counting and/or quantitative determination of hemoglobin are carried out directly

## Principle

It lyses red blood cells while keeping the morphology of white blood cells unchanged, enabling the analyzer to count and classify white blood cells by the laser light scatter method...

