


In Vitro Diagnostic Reagent

Diluent-UT II

Used to dilute the sample and prepare a cell suspension prior to blood cell analysis

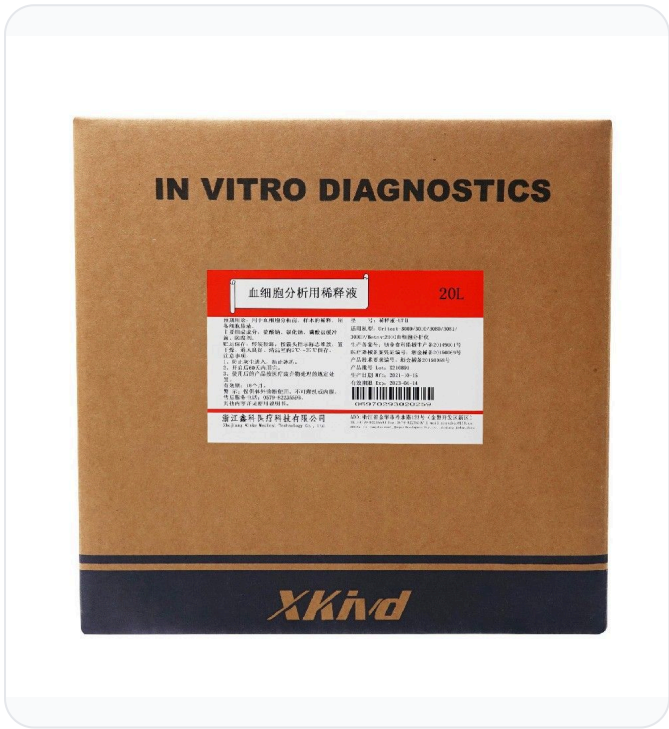

Product Photo

Product Model	Diluent-UT II
Packaging Specifications	5 L, 20 L
Product Code	p1191

Intended Use

Used to dilute the sample and prepare a cell suspension prior to blood cell analysis

Principle

It dilutes the blood cells and, in combination with the flow cytometric technique based on the impedance method, colorimetry, laser light scatter and fluorescence staining, performs cell differentiation, counting and hemoglobin determination

Composition

Sodium sulfate, sodium chloride, phosphate buffer solution, preservatives...

This document is compiled from product registration information for reference only. Refer to the product instructions and the registration certificate for details.